smartLAB[®] scale

Personenwaage

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Bathroom Scale

User Manual

Please read the user manual carefully before the first use.

Current information about this device, visit smartLAB on

www.support.hmm.info

HMM Diagnostics GmbH
Friedrichstraße 89
D-69221 Dossenheim, Germany
www.hmm.info

Content

I. Introduction	20
Principle of Measurement	20
II. Your smartLAB® scale W	21
Display & functions	21
Content	22
General Information	22
Security Information	23
III. Setup & Operating functions	24
Inserting and Changing the batteries	24
Initial start - „quick weighing “	25
Start Measurement	25
Optional: Weighing on a carpet	26
IV. Miscellaneous	27
Maintenance	27
Troubleshooting	27
Symbol Meaning	29
Regularity and Safety Notice	30
Warranty	30

I. Introduction

Thank you for choosing our **smartLAB**[®] *scale* Bathroom Scale. This product enables you to measure your body weight. Please read the instructions carefully before you first use the scale and keep the instructions for a future use.

Principle of measurement


Your **smartLAB**[®] *scale* Bathroom Scale measures your body weight with the help of 4 high accurate sensores.

II. Your smartLAB[®] scale

General Instructions

The **smartLAB[®] scale** Bathroom Scale applies SENSE ON technology. It's not necessary to turn on or turn off the scale, just stand on the scale for measurement. It will turn on automatically when stepping on.

Display & functions

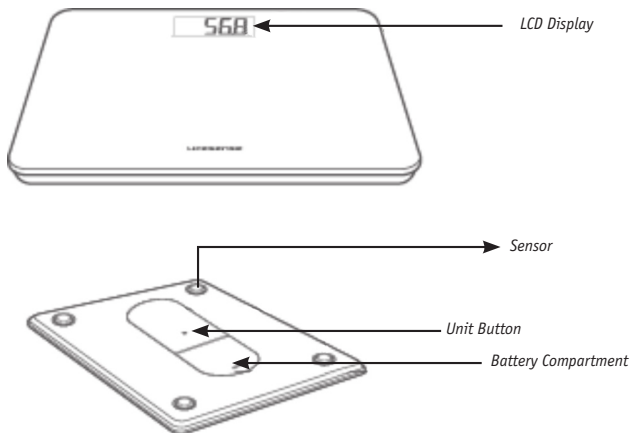LCD Display

kg	Kilogram
lb	Pound
	Low Battery

Specifications

1. Device: **smartLAB**[®] scale
2. Max. capacity: 150 kg / 330 lb
3. Units: kg / lb (in 100 g / 0,2 lb steps)
4. Buttons: Units
5. Battery: 4 x AAA batteries (1,5 V)
6. Size (L x W x H): 321 x 321 x 23,5 mm
7. Weight: 1980 g (without batteries)
8. Operating temperature: 0°C - 40°C
9. Storage temperature: -20°C - 60°C
10. Symbol when overburdened: „----“
11. Symbol for low batteries: „“
12. Communication/Data transfer: ANT+ and BLE wireless

Contents

- 1 **smartLAB**[®] scale Bathroom Scale
- 1 User manual
- 4 1,5V „AAA“ batteries
- 4 carpet extension feet

General Information

- Please place the scale on an even, flat and hard floor. The floor is the basis for a correct measurement. Soft surface such as carpet will affect the performance of the scale. Measurement on a carpet can be enabled by the extension feet.
- If possible please use the scale at the same time every day (we recommend in the morning), after having used the toilet, with empty stomach and without clothes in order to be able to compare the results.
- Please wait about 15 minutes after getting up before you measure your body weight so that the water of your body can disperse.
- Please do not move during the measurement.
- Please wait several hours before measuring after unaccustomed physical exercise, sauna or bath, drinking and eating.

- The device is for your own use only and not for medical or commercial use.
- Please note that technical related measuring tolerances are possible since the scale is not especially made for professional medical use.
- The **smartLAB**[®] *scale* Bathroom Scale corresponds to the EC directive 89/336.

Please feel free to contact the manufacturer or the customer service for more information.

Security information

-  Avoid slippery floor such as tile floor
-  Avoid using the device close to a microwave oven
-  Please do not use the **smartLAB**[®] *scale* Bathroom Scale with wet feet – slip hazard! Avoid measuring platform immediately after bath.
-  Do not swallow the batteries – danger of life. Please store batteries and Bathroom Scale at a place children cannot reach. In case a battery has been swallowed, please consult a physician immediately.
-  Keep children away from the packaging materials - danger of asphyxiation!
-  Do not recharge the batteries. Do not take them to pieces. Do not throw them into a fire. Do not short-circuit them.
-  Avoid storage in following locations:
 - Where there is water
 - Where the device may be exposed to extreme temperature, humidity, moisture, direct sunlight, dust or salt air
 - Where there is risk of shock or drop
 - Where you store chemicals or full of corrosive gasses

III. Setup & Operating functions

Inserting and changing the batteries

- Please open the cover of the battery case at the back of the device.
- Insert 4 AAA batteries (1,5 V). Please pay attention to the polarity.
- The digits “8888” will appear on the LCD Display

- Close the cover of the battery case and wait until the digits “0.0kg” show on the display.

Note: In case the scale does not operate please remove the batteries and insert them again.

 Your **smartLAB**[®] *scale* Bathroom Scale has a symbol indicating that the batteries need to be changed. In case the batteries are almost empty, the symbol “” will show on the display. The device will power off in 4 seconds. You should change the batteries then.

 Please remove the batteries in case you are not going to use the scale for one month or more in order to avoid an automatic discharge of the batteries.

The Bathroom Scale as well as the batteries must be disposed according to the local regulations.

 In order to guarantee a correct measurement you should calibrate the scale always after any location change. Just tip on the scale and wait until it shows “0.0” on the scale.


Initial Start - „quick weighing“

Use the **smartLAB®** scale Bathroom Scale to measure your weight. Do not move, stand still and stand on both legs. The scale will turn on automatically when stepping on it. The scale starts measuring immediately, a beep tells you the end of the measurement and finally shows your body weight on the display. Now you can step off the scale. The scale will deactivate automatically after several seconds.

Start Measuring

Select Measuring Unit

With batteries correctly installed, press the “UNIT” button on the back of the scale to select the measurement unit. The default measurement unit is “kg”. You may press “UNIT” button to choose among kilogram and pound.

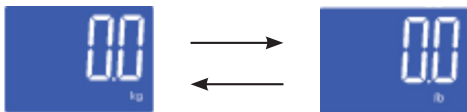Daily Measurement

Thanks to the SENSE ON patent technology, the scale will automatically turn on when you step on the platform.



Stand still and keep full contact with the scale until the LCD displays the weight and the unit stops blinking. A beep will let you know when the weighing is done. The scale will turn off automatically after a few seconds.

Optional: Weighing on a carpet

You can use the **smartLAB**[®] scale with the extension feet on a carpet. For this please follow these steps:

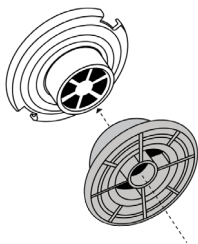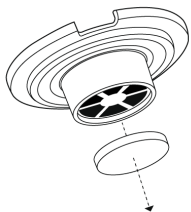1. Remove the non-slip pads from the scales feet
2. Fit the extension feet firmly onto the scales feet. Now you can use the scale on a carpet.

IV. Miscellaneous

Maintenance

- Please clean your device regularly with a wet soft cloth. Use a dry soft cloth to dry up the device. Do not use caustic cleaning solutions. Do not put the device under fluent water.
- Protect the device from impacts, humidity, dust, chemicals or fluctuations of temperature. Do not store near heat sources such as ovens or heating installations.
- NEVER take the scale to pieces by yourself. In this case the warranty will not be valid anymore. Only the manufacturer or the authorized trader is allowed to repair the scale. Please check the batteries before any complaint and change them if necessary.

Troubleshooting

Error Prompt

Error	Description	Solution
	Overload. The device will power off in four seconds.	Stop using this scale for measurement. The scale only weighs until 150kg
	Low Battery. The device will power off in four seconds.	Replace all four batteries at the same time. Please purchase the authorized batteries for replacement.

Problem	Root Cause	Solution
Abnormal measuring results: - Too high, or - Too low, or - Huge difference between two recent measurement	Incorrect posture	Please step on the platform and stand still.
	The device is located on a soft ground such as a carpet OR on a rugged surface.	Please place the device on a flat, hard surface.
No display on LCD when the device powers on	Batteries not yet installed.	Install the batteries. (Please refer to "Insert the batteries")
	Worn batteries.	Replace all four batteries at the same time. Please purchase the authorized batteries for replacement.
The device powers off automatically.	Low battery.	Replace all four batteries at the same time. Please purchase the authorized batteries for replacement.

Symbol Meaning

The warning signs and symbols are essential to ensure your correct and safe use of this product and protect you and others from injury. Please kindly find the meanings of the warning signs and symbols, which you may encounter in the label and user manual, as follows:

	Symbol for "THE OPERATION GUIDE MUST BE READ"		Symbol for "MANUFACTURER"
	Symbol for "COMPLIES WITH MDD93/42/ECC REQUIREMENTS"		Symbol for "MANUFACTURE DATE"
SN	Symbol for "SERIAL NUMBER"		Symbol for "ENVIRONMENT PROTECTION – Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice"
	Symbol for "DIRECT CURRENT"		

WEEE note

The WEEE (Waste Electrical and Electronic Equipment) Directive, that came into effect as European law on 13th February 2003, led to a major change in the disposal of electrical equipment.

The primarily purpose of this directive is to prevent electronic waste (WEEE), while encouraging for reuse, recycling and other forms of reprocessing to reduce waste.

The WEEE logo on the product and the packaging indicates that it is not allowed to dispose of the product in the ordinary household waste. It is your responsibility to deliver all the disused electrical and electronic equipment to the respective collection points. A separate collection and proper recycling of electronic waste helps dealing economically with natural resources. Furthermore, the recycling of electrical waste is a contribution to the conservation of our environment and thus the health of all people. More information concerning the disposal of electrical and electronic devices, reprocessing and the collection points

you can get from the local authorities, waste disposal companies, from retailers and manufacturers of the device.

RoHS compliance

This product complies with Directive 2002/95/EC of the European Parliament and the Council of 27th January 2003 regarding the limited use of dangerous substances in electrical and electronic equipment (RoHS) and its variations.

Warranty

HMM Diagnostics GmbH products need to fulfill high quality requirements. Because of this reason, HMM Diagnostics GmbH gives a 2-year warranty by purchasing this **smartLAB**[®] product. You can even extend the warranty from 3 to 5 years without extra pay when you register your product. Please register online under following website: www.hmm.info/en/registration.

Wear parts, batteries etc. are excluded from warranty

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